


SIRABE

MICRO CATHETER SYSTEM


PURE SELECTIVE

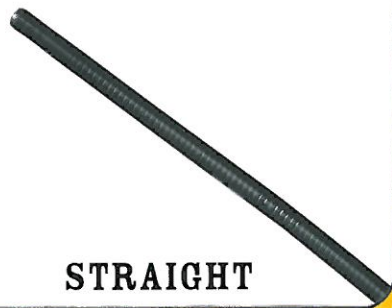

STRAIGHT



SPIRAL

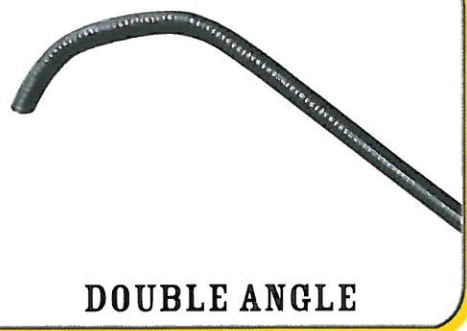


DOUBLE ANGLE

PIOLAX

Product characteristics

- **“Silicone-coated PTFE liner”**
Easier insertion through the shaft lumen
 - Smoother insertion of embolization coils, and improved operability and trackability of microguidewires.
- **Unique tip shapes: straight, spiral and double-angled**
 - 3D-spiral and 2D-double-angled tips can accommodate difficult cases.
 - Tip shape can be modified using steam (with 90° mandrel angle, shaping angle is about 60°)
- **“Velvet Feel”**
Proper shaft balance and anti-kink tip
 - The “anti-kink soft tip” design at the tip has good trackability.
 - The flexible tapered shaft needs less back-up force of guiding catheter and is suitable for use with highly tortuous vessels.
 - The rigid handle shaft enables safe manipulation and allows superior pressure application.
- **Highly selective performance**
2.2Fr. narrow shaft and hydrophilic surface coating
 - Ensures smooth access to more distal vessels.
 - High-durability coating achieves stable lubrication, even with frequent usage.

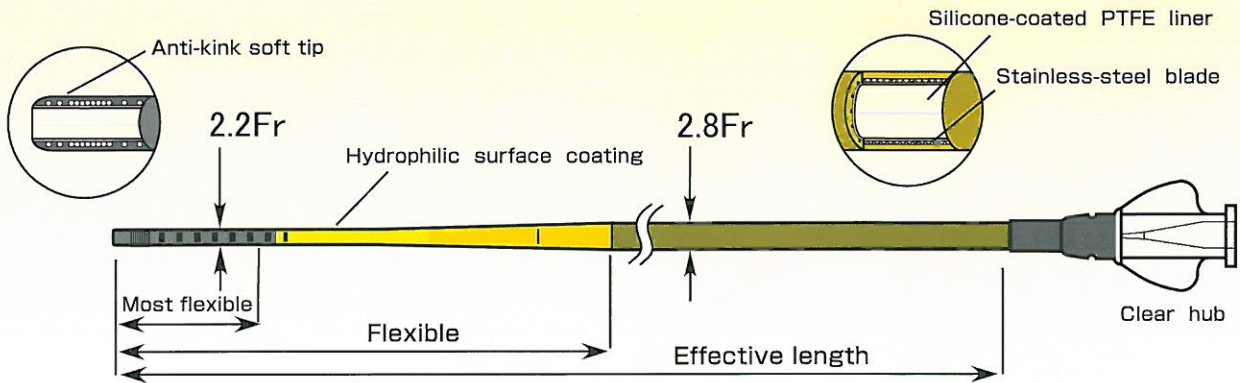

Catarg No.	Tip shape	External diameter (tip/handle)	Internal diameter	Effective length/Flexible length/ Most-flexible length	Guaranteed compressive strength	Hydrophilic coating length	JAN code
MC18-S110K	Straight	2.2/2.8Fr (0.74/0.94mm)	0.021 inch (0.54mm)	110cm/38cm/25mm	800psi (5520kPa) ※with open top	65cm	4545428010485
MC18-S130K				130cm/38cm/25mm			4545428010492
MC18-W110K	Spiral			110cm/38cm/25mm			4545428010546
MC18-W130K				130cm/38cm/25mm			4545428010553
MC18-D110K	Double-angled			110cm/38cm/25mm			4545428010515
MC18-D130K				130cm/38cm/25mm			4545428010522

- Approval name: マイクロカテーテル
- Approval number: 21700BZZ00138000
- Sterilized with ethylene oxide ● Do not reuse ● Medical device classification: Highly regulated equipment



Manufacturer :
Piolax Medical Devices, Inc.
179 Kariba-cho Hodogaya-ku Yokohama,
240-0025, Japan
TEL.+81-45-710-1570 FAX.+81-45-710-1584